nanosonics
Infection Prevention. For Life.

trophon® Printer
User Manual

ENGLISH **2**
DEUTSCH **19**
FRANÇAIS **36**
NEDERLANDS **53**
ITALIANO **70**
NORSK **87**
SUOMI **104**
DANSK **121**
SVENSKA **138**
ESPAÑOL **155**
TÜRKÇE **172**
日本語 **189**

ENGLISH

nanosonics
Infection Prevention. For Life.

trophon® Printer
User Manual

Read this manual before operating the trophon® Printer to ensure correct procedures are followed.

For further information, contact your customer service representative or visit the Nanosonics website.

All technical specifications and system approvals are listed in the section: Printer Technical Specifications N05036/DTP200.

© 2026 Nanosonics Limited. All rights reserved.

The contents of this manual are correct at the time of product sale.

Nanosonics' technologies are protected by a combination of patents, trademarks, and exclusive proprietary rights.

trophon® is a registered trademark of Nanosonics Limited.

nanosonics
Infection Prevention. For Life.

www.nanosonics.com

CE Conforming to Low Voltage Electrical Equipment Directive (EU) 2014/35/EU & Electromagnetic Compatibility Directive 2014/30/EU

UKCA Conforming to UK Electrical Equipment Safety Regulations 2016 (SI 2016/1101) and UK Electromagnetic Compatibility Regulations 2016 (SI 2016/1091)

Table of Contents

ENGLISH 2

INTRODUCTION 5

- **Labels, Symbols and Warning** 5
- **What's in the Box** 6
- **Printer Features** 7

INSTRUCTIONS FOR USE 9

- **First Use and Setup** 9
- **Installing the Label Roll** 10
- **Confirming Printer Operation** 12

ROUTINE CARE AND MAINTENANCE 13

- **What You Will Need** 13
- **Warnings** 13
- **How to Clean the Printer** 13
- **Cleaning Schedule** 15

TROUBLESHOOTING 16

- **Connecting and Disconnecting the Printer** 16
- **Printer – Error Messages** 16
- **Troubleshooting Steps** 16
- **Printer Disposal** 17
- **Printer Technical Specifications N05036/DTP200** 17

INTRODUCTION

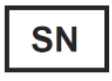
The trophon Printer provides a label record of the most recent high level disinfection cycle directly from the trophon device. The label is generated at the completion of each high-level disinfection cycle and can be placed directly on to a patient medical record or other documentation. Printed label records remain legible for up to 7 years if stored as per the specifications outlined on the trophon printer label roll packaging.

This printer is compatible with the trophon2 and trophon3 devices only. Accuracy and lifetime of records is not guaranteed if an unapproved accessory (printer or label roll) is used or if records are stored outside the storage specifications.

For any operation, fault or maintenance queries, please contact your customer service representative.

Labels, Symbols and Warning

The following warnings and symbols apply to the trophon Printer

	Caution		Consult Instructions For Use
	Handle With Care		Electrical Shock Hazard
	Cut Hazard		Moving parts can crush/cut
	Fragile		This Side Up
	Keep Out of Direct Sunlight		Avoid Moisture
	Legal Manufacturer		Date of Manufacture
	Warning: Hot Surface		Warning: Moving parts
	Separate collection for electrical and electronic equipment		Recyclable Corrugated Cardboard
	Conforming to Low Voltage Electrical Equipment Directive (EU) 2014/35/EU & Electromagnetic Compatibility Directive 2014/30/EU		RoHS Compliant (Restriction of Hazardous Substances)
	Regulatory Compliance Mark		
	Registration, Evaluation, Authorisation and Restriction of Chemicals European Union regulation dating from 18 December 2006, amended on 16 December 2008 by Regulation (EC) No 1272/2008.		
	Product Number		Serial Number

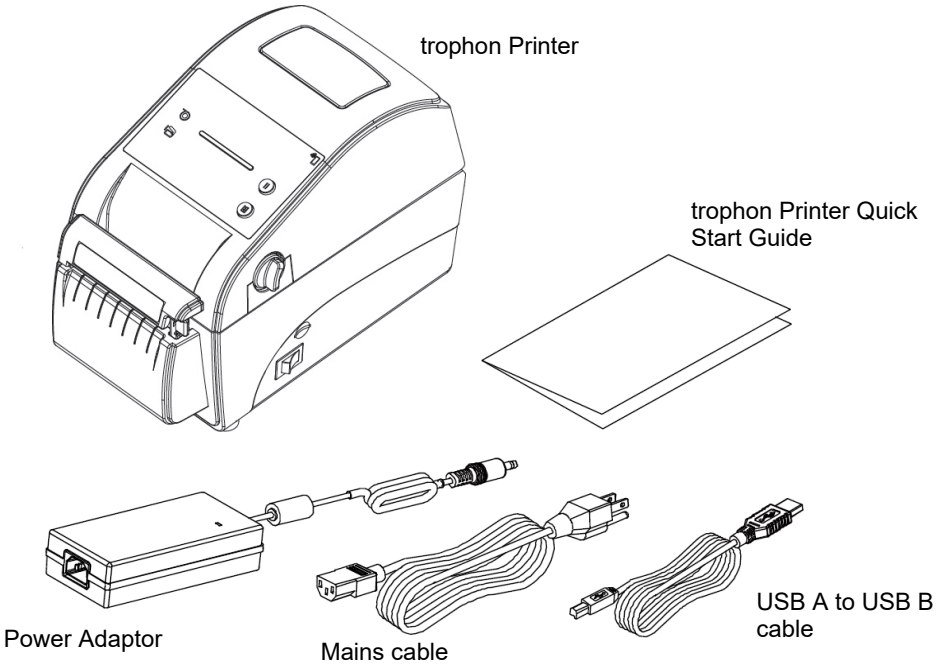
UK CA	Conforming to UK Electrical Equipment Safety Regulations 2016 (SI 2016/1101) and UK Electromagnetic Compatibility Regulations 2016 (SI 2016/1091)
TÜV SUD C US	UL 62368-1:2019 CSA C22.2 No. 62368-1:2019
FC	Complies with FCC Part 15, Subpart B (Class B device)
FSC MIX Packaging \| Supporting responsible forestry FSC™ C156753	Forest Stewardship Council for sustainable forestry FSC Mix can be a blend of FSC-certified, recycled and controlled sources
energy ENERGY STAR	Energy Star

What's in the Box

The trophon Printer has been packaged to withstand damage during shipment. Where possible, keep the packaging materials in case the printer needs to be returned. When unpacking the printer ensure that you have received the items shown in Figure 1 .

NOTE: If any parts are missing or damaged, contact your customer service representative.

Figure 1 – trophon Printer Box Contents

Printer Features

The main components and their locations on the trophon Printer are shown below in Figure 2 to Figure 5 .

Figure 2 – trophon Printer External Features

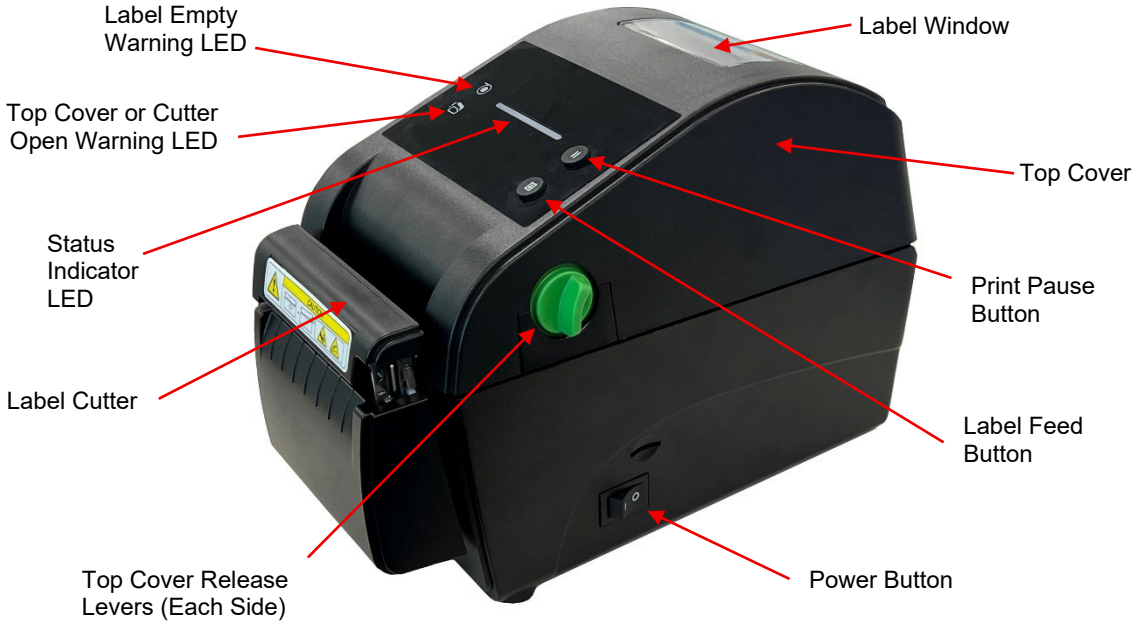
Figure 3 – trophon Printer Internal Features and Connectivity

Figure 4 – trophon Printer Control Panel

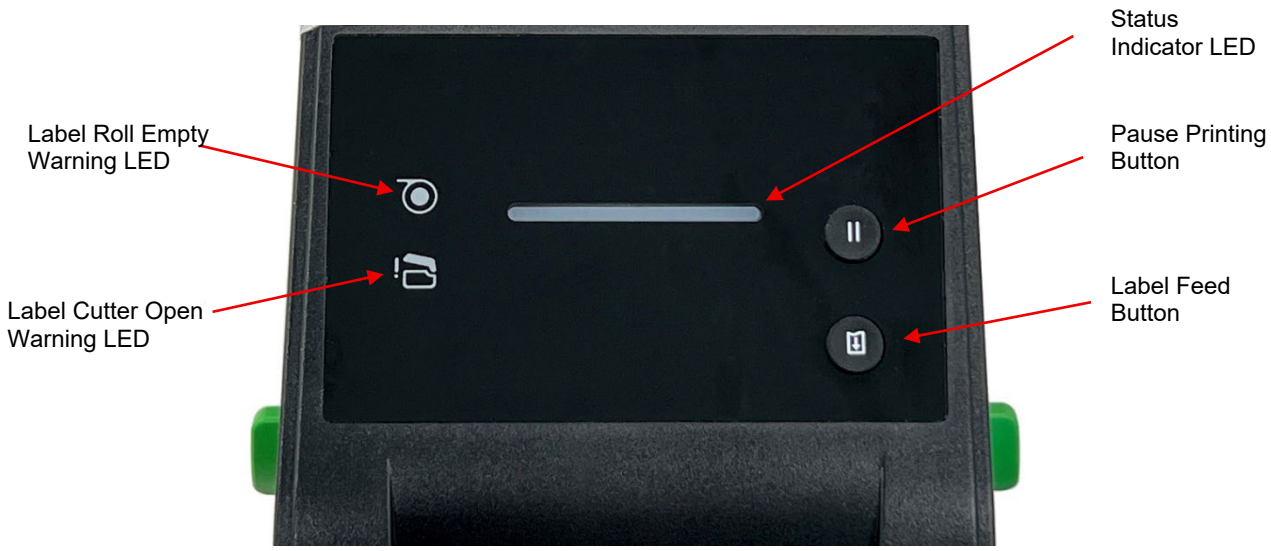
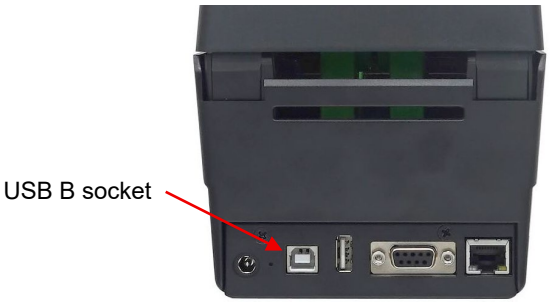
Figure 5 – trophon Printer Status Indicator LED Details

Green	**ON:** Printer is ready **Blinking:** Printer job downloading or printer is paused
Amber	The system is busy
Red	**ON:** Printer top cover or label cutter is open **Blinking:** System error: label roll jammed or empty

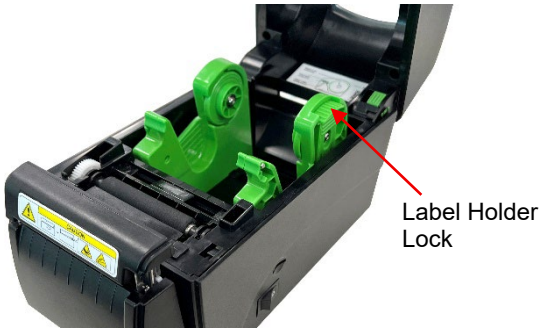
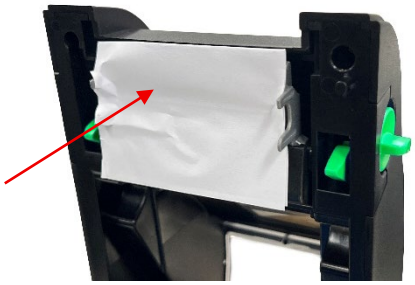
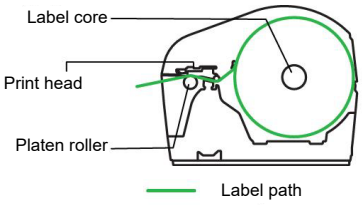
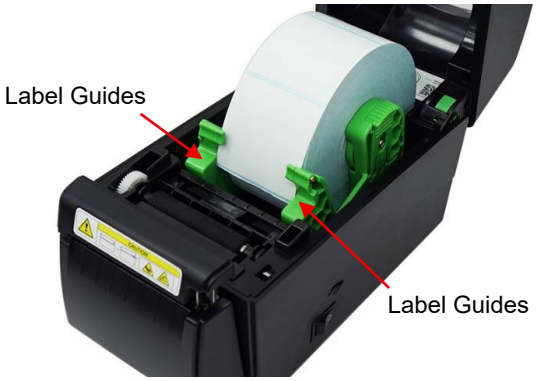
INSTRUCTIONS FOR USE

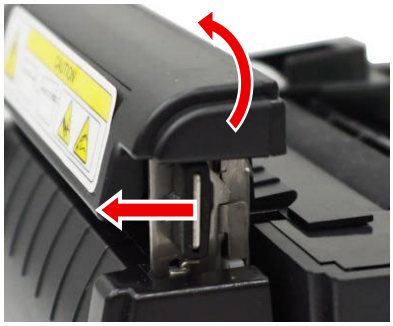
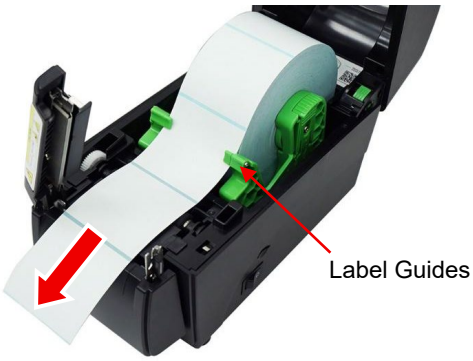
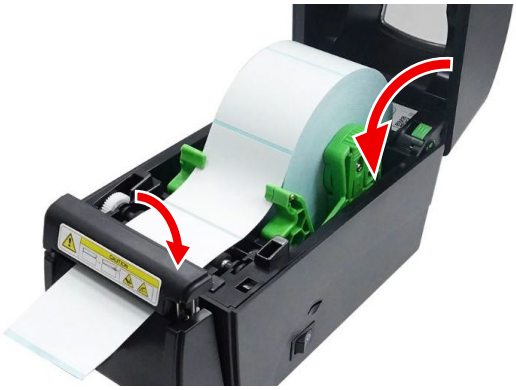
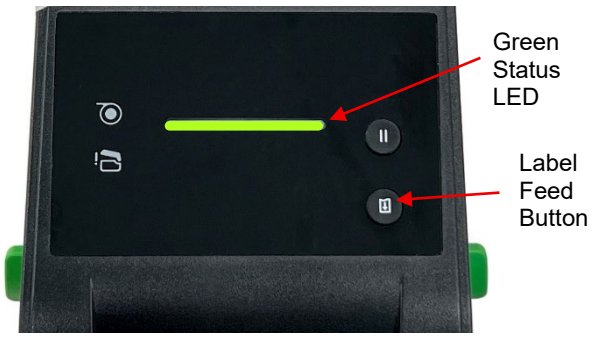
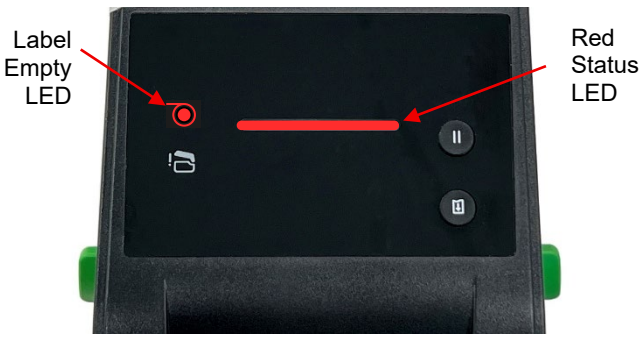
First Use and Setup

1. Unpack the printer from its packaging and ensure all the parts are present, see Figure 1 . Remove all adhesive tape and protective materials.	
2. Position the printer securely on a bench beside the trophon device. The trophon Printer Wall Mount or trophon Printer Cart Mount accessories can also be used if the printer needs to be fixed to a wall or used on the trophon cart.	
3. Ensure the power switch on the right side of the printer is set to OFF.	
4. Connect the power cord to the power adapter.	

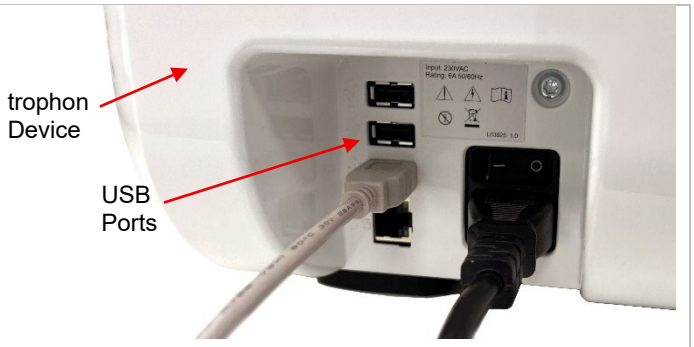
5. Connect the power adapter to the power jack at the rear of the trophon Printer. Fully insert the power cord plug into the power outlet socket and turn it on.	
6. Connect the printer to your trophon using the supplied USB A to USB B cable.	

Installing the Label Roll

1. Open the top cover and pull apart the label holders. Press the label holder lock down to lock the label holders open.	
2. Remove the protective paper.	
3. Insert the label roll with the printable side facing upwards and release the label holder lock to hold the label roll. When installed, ensure the label roll can rotate freely. 	

4. Open the label cutter gate by pulling the release button towards the front of the printer to unlock it.	
5. Thread the label under the label guides and over the platen roller as shown in the diagram below.	
6. Close the label cutter and make sure it clicks into position. Close the top cover.	
7. Turn the printer on and wait for the green status LED to light. Press the label feed button to check the label is loaded correctly and feeds properly. The label cutter should also cut the labels but leave a small piece allowing the label to be torn off by hand.	
8. If the printer does not turn on or the status LED is red, refer to the Troubleshooting Section During normal operation the trophon device will display an onscreen notification when the label roll requires replacement, and the printer will display an error, and the label empty warning LED will be lit.	

Confirming Printer Operation

1. Once the printer is connected to the trophon device via USB cable and is powered on, the printer will be automatically detected by the trophon device.	
2. In the trophon user interface navigate to: **Main Menu** > **Settings** > **Printer Label** > **Test Print** OR **Main Menu** > **Shortcuts** > **Printer Label** > **Test Print**. Select the test print option and ensure the trophon Printer is working correctly. If not, refer to the Troubleshooting section.	

ROUTINE CARE AND MAINTENANCE

What You Will Need

The following supplies are recommended for cleaning the printer:

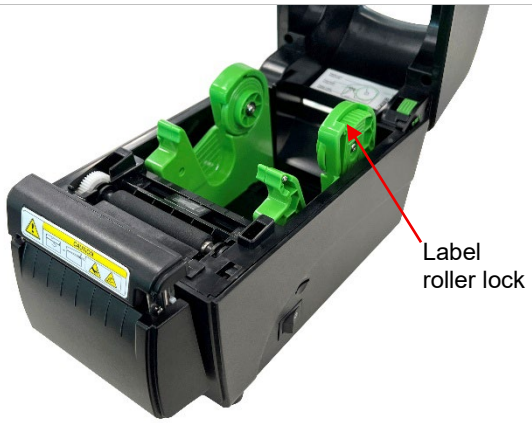
- Cotton swabs
- Quaternary ammonium (quat) wipes for disinfecting the printer exterior
- Soft brush with non-metallic bristles
- Isopropyl alcohol

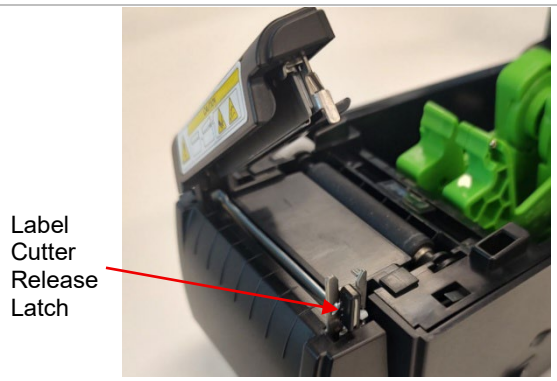
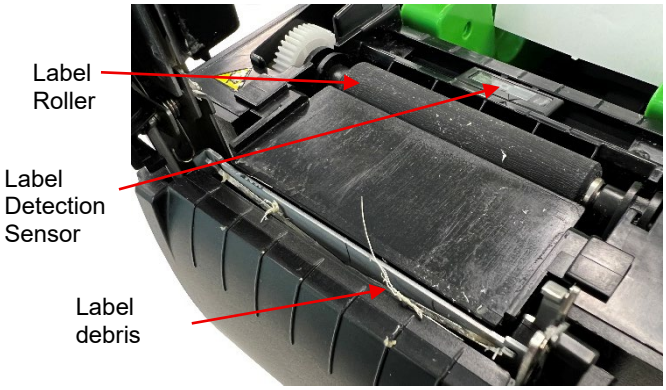
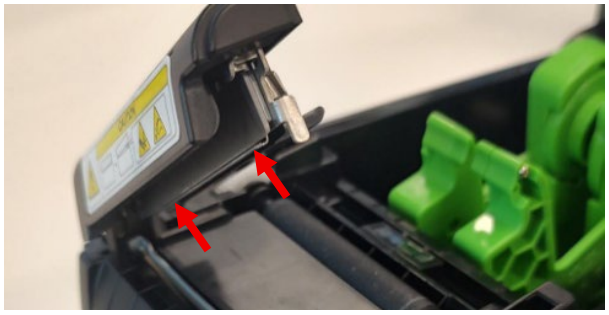
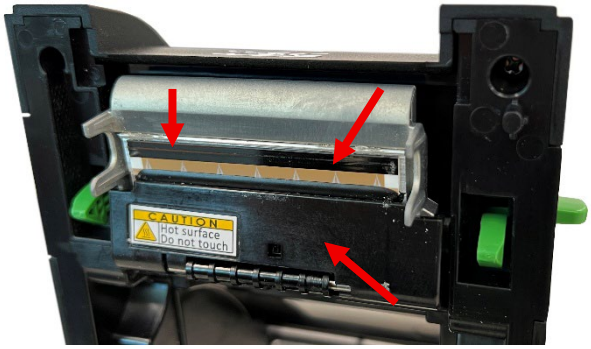
Warnings

The inside of the printer may collect label dust, adhesive residue etc during normal operation. To ensure print quality and reliability the printer should be cleaned periodically as shown in the Cleaning Schedule.

- The printer can be cleaned and/or disinfected if necessary, using isopropyl alcohol or quaternary ammonium (quat) wipes
- Only use the cleaning materials recommended in this document. Use of other materials may damage the printer and void its warranty.
- Do not spray or drip liquid cleaning solutions directly into the printer. Apply the solution to a cotton swab or use a cleaning wipe.
- Do NOT submerge the printer or pour liquids over it. Keep the power jack completely dry.
- Do not use pressurised (canned or compressed) air to clean the interior of the printer as it can blow dust and debris onto sensors and other critical components.

How to Clean the Printer

1. Turn the trophon Printer off but leave it connected to the power adapter. If the printer has been used, allow the printhead to cool for 1 minute.	
2. Remove any metal objects from your hands, i.e. rings or bracelets. Open the top cover and remove the label roll. Lock the label roller open using the label roller lock.	

3. Open the printer top cover and open the label cutter by pulling the label cutter release latch.	
4. Check the label roller, label cutter and sensor for debris. Carefully remove any large pieces of label from the cutter using tweezers or similar tools. Use a cotton swab containing isopropyl alcohol to clean these components and remove any adhesive residue.	
5. Clean the cutter blade with a cotton swab containing isopropyl alcohol.	
6. Clean the printhead when replacing the label roll using a cotton swab containing isopropyl alcohol. Do not touch the printhead except with a cotton swab.	
7. For external cleaning and disinfection, wipe the outside covers of the printer with isopropanol or quaternary ammonium (quat) wipes until all surfaces are visibly clean. Ensure that liquid does not contact the power socket, printer interior, or label roll.	

Cleaning Schedule

Component	Recommended Schedule
Printhead	Clean the printhead when you load new labels
Platen Roller	Clean the platen roller when you load new labels
Label Cutter	Clean as needed
Label Sensor	Clean the sensor monthly
Exterior	Clean as needed
Interior	Clean as needed

TROUBLESHOOTING

Connecting and Disconnecting the Printer

The printer is automatically detected and enabled when connected to the trophon device.

The printer can be enabled or disabled from the printer settings accessed on the trophon device screen after it has been plugged into/unplugged from the trophon device.

Printer – Error Messages

Printer error messages are displayed on the screen of the trophon device. If an error message is displayed, two options are available:

Retry

- After the error has been addressed, select Retry. The trophon device will attempt to re-print the last disinfection record.
- If the last disinfection record did not reprint, the Retry and Ignore options will be offered again. You can attempt to fix the printer issue and then attempt reprinting.

Ignore

- If the printed record is not required and you do NOT want to *Retry* the print, select *Ignore.* If Ignore is selected, the disinfection label record will not be printed.

Troubleshooting Steps

Problem	Possible Cause	Recommended Solution
Printer does not turn on. Power indicator is not lit	• The power cord is not properly connected.	• Check the power cord is securely plugged into the printer and outlet. • Ensure printer is switched on
Top Cover Open LED is lit	• The top cover or label cutter is open.	• Close the top cover and label cutter are both closed and clicked into place.
Not Printing	• USB cable is not connected properly or is damaged.	• Ensure both ends of the USB A to USB B cable are connected to the back of the trophon Device and the trophon Printer • If the problem continues, try a different USB A to USB B cable.
No print on the label	• Label roll is not loaded correctly or is the wrong type. • Printhead is dirty	• Ensure the label roll is loaded correctly and genuine Nanosonics trophon Printer Label Rolls are used • Clean the printhead.
Label not feeding	• The label is installed incorrectly.	• Check the operation of the label roll. • Install a new label roll.
Label jam	• Ensure label roll is loaded correctly. • Labels may be stuck inside the printer mechanism.	• Ensure the label roll is threaded correctly and passes under both label guides. • Remove any labels stuck inside the printer mechanism.
Poor print quality	• Label roll is loaded incorrectly • There is dust or adhesive on the printhead. • Printhead element is damaged.	• Check and/or reload the label roll. • Clean the printhead and platen roller. • The label release lever is not latched properly.
Gray line on the blank label	• The printhead is dirty. • The platen roller is dirty.	• Clean the printhead. • Clean the platen roller.
Label feeding is not stable (skewed) when printing	• The label guides do not touch the edge of the label roll.	• Open the top cover and ensure the label is passing under the label guides on both sides
Damaged Label Roll	• Label roll installation or feeding is incorrect. • Incorrect label type.	• Make sure the label roll touches the edge of the label guides. • Ensure you only use trophon Printer Labels.

Printer Disposal

The following information is for EU-member states only:

This symbol indicates that the product must not be treated as general waste and must be disposed of correctly. For more detailed information about the collection and recycling of this equipment, please contact your supplier.

REACH COMPLIANCE STATEMENT

The Registration, Evaluation, Authorization and Restriction of Chemicals Regulation EC 1907/2006, commonly referred to as REACH, is a European regulation on chemicals and their safe use. With the publication of candidate lists for authorization, the European Chemicals Agency identifies sets of substances of very high concern (SVHC) that manufacturers must disclose to customers if used in articles in their products above 0.1% by weight. Based on the current version of the candidate list for authorization (250 chemicals), this product contains the following SVHC above the disclosure threshold:

SVHC	CAS Number	Application
2,2',6,6'-tetrabromo-4,4'-isopropylidenediphenol	79-94-7	ABS plastic

The SVHCs used in this product does not pose a safety risk to customers under normal use conditions

Printer Technical Specifications N05036/DTP200

N05036 Electrical Specification –printer	Input: DC 24V, 2.5A
Electrical Specification – External universal switching power supply	Input Voltage: 100 – 240Vac Input Current: 1.8A, 50/60Hz Mains Inlet: IEC type C13 Equipment must be connected to an earthed outlet using power cable supplied. Output: DC 24V, 2.50A
Environmental Specification	Operation: 41°F to 104°F (5°C to 40°C), 25 to 85% non-condensing Storage: -40°F to +140°F (-40°C to +60°C), 10 to 90% non-condensing
Physical Dimension	138 mm (width) × 172 mm (height) × 210 mm (depth) 5.43” (width) × 6.77” (height) × 8.27” (depth)
Data Interface	USB B port
Electromagnetic Compliance	FCC Class B, CE, AS/NZS CISPR 32:2015, ISED (ICES-003, Issue 7, October 15, 2020 Class B) This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures: — Reorient or relocate the receiving antenna. — Increase the separation between the equipment and receiver. — Connect the equipment into an outlet on a circuit different from that to which the receiver is connected. — Consult the dealer or an experienced radio/TV technician for help.

	This equipment complies with Part 15 of FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. This Class B digital apparatus complies with Canadian ICES-003. Cet appareil numérique de la classe B est conform à la norme NMB-003 du Canada.

DEUTSCH

nanosonics
Infection Prevention. For Life.

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Nanosonics, Inc.
Nanosonics Europe GmbH
Nanosonics UK Limited